

Vaporizer

Module Monitoring




Maintenance-free Vaporizer

● Vaporizer

● Mindray V60 vaporizer (Iso, Sev)

- ✓ Automated flow/ Temp/ Pressure compensation ensure a constant output and accurate concentration under different environment conditions.
- ✓ Good air-tightness, reduce agent gas waste
- ✓ **Maintenance-free** design ensure that no calibration is required throughout the entire life span
- ✓ **FDA Approval**
- ✓ Large capacity (300 ml with most wick, 360 ml with dray wick)
- ✓ Highly compatible: support three different filling options and standard selectatec mounting system.

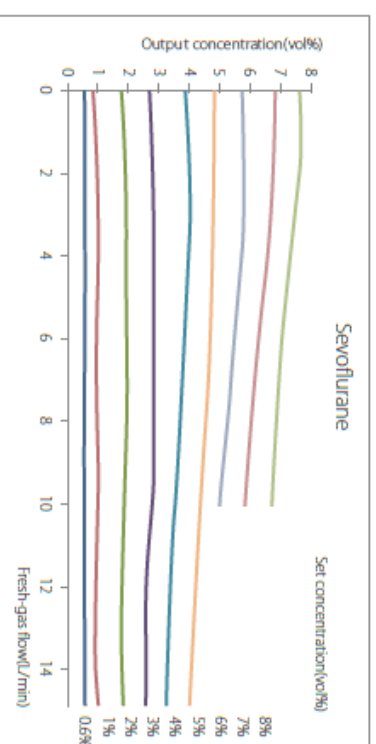

The diagrams show the influence of flow on the concentration delivered after 1 minute at 22°C, 101.3 kPa when operating with Air.

Vaporizer

Pic	Part number	Description	Apply to
	115-005348-00	Mindray V60 vaporizer, Isoflurane, Pour-fill	WATO
	115-005345-00	Mindray V60 vaporizer, Isoflurane, Key-fill, with Key-fill adaptor	
	115-005349-00	Mindray V60 vaporizer, Sevoflurane, Pour-fill	WATO
	115-005346-00	Mindray V60 vaporizer, Sevoflurane, Key-fill, with Key-fill adaptor	
	115-005350-00	Mindray V60 vaporizer, Sevoflurane, Quik-fill	

- This part is for WATO series and fits Selectatec manifold.
- At present, the Veta vaporizer is fixed in the main unit, and is the service part. please contact the service team to order.



Monitoring Sensor

Pic	Part number	Description	Apply to
	115-001366-00	Flow sensor kit, including: -0601-30-69700 Flow sensor, inspiration, 1 pc -0601-30-78894 Flow sensor, expiration, 1 pc	WATO

● FlowSecure Patent Technology, What you set is what you get

- ✓ Dynamic tidal volume compensation, prevent the influence of fresh gas flow or system compliance & leakage
- ✓ High-quality insp & exp flow sensors, precise and long service life
- ✓ Internal the 3rd flow sensor for calibration and safety backup

